
truelife

OXIMETER X5 BT

Pulse oximeter

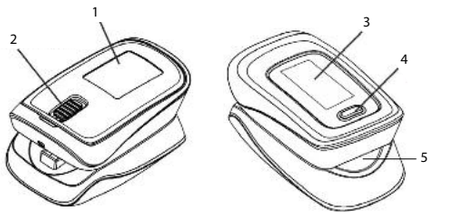EN – USER MANUAL

1. BOX CONTENTS

1× TrueLife Oximeter X5 BT
2× AAA alkaline battery
1× Suspension cable
1× Case
1× Manual
1× Safety instructions


The TrueLife Oximeter X5 BT uses a light source that passes through the tissue and hits a detector that measures its intensity. The Pulse Oximeter X5 BT measures oxygen saturation or blood oxygen levels as well as heart rate.

2. PRODUCT DESCRIPTION


1. Battery cover
2. Battery cover lock
3. Display
4. Control button
5. Finger clamp

Display Description

1. Oxygen saturation symbol
2. Oxygen saturation (percentage value)
3. Pulse rate waveform
4. Battery status symbol / Bluetooth symbol (alternately)
5. Heart rate symbol
6. Perfusion index
7. Pulse rate
8. Pulse rate bar

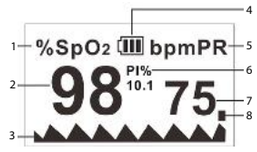
3. BATTERY INSTALLATION

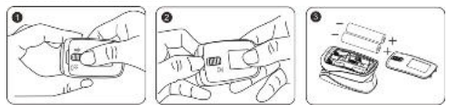
- Remove the battery cover on the back of the oximeter and insert 2 AAA batteries. Make sure that the batteries are inserted with the correct polarity as indicated in the battery compartment. Do not use batteries after their expiration date.
- If the battery status symbol flashes on the display, the batteries must be replaced. When battery power is critically low, the device will not take a reading.
- If you will not be using the device for an extended period of time, remove the batteries.

4. SETTINGS

- The oximeter offers the possibility to set audible a alert, pulse rate beep and brightness. You can then set the upper and lower limits for the audible alert.
- To enter settings, turn on the oximeter by pressing the control button. After the display turns on, hold down the control button again for 2 seconds. The display shows the menu for setting the alert, the pulse rate beep and the brightness. The menu for setting the alert limits is accessed by holding down the control button while the asterisk is next to „Ale setup".

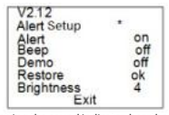Setting the sound indicator, heartbeat soundtrack and brightness

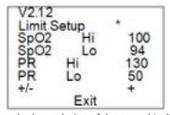Setting the boundaries of the sound indicator

- To exit settings, navigate to the line „Exit" and hold down the control button .

Setting the measurement alert

- The reading can be accompanied by an audible alert.
- To turn on the alert, scroll to the „Ale" line and hold down the control button until „on" is set. To turn off the alert, proceed in the same way and set to „off".

Setting audible alert limits

- When taking a reading, the oximeter starts emitting an audible alert when one of the measured values is smaller or larger than the set limits.
- To change the limits, enter the alert limit menu and select the line with the value you want to change. Once you are at your desired value, hold down the control button and the value will change. If „+" is on the bottom line, the value will rise. If „-" is there, the value will decrease.

Note: We recommend that you leave the alert limits at the preset values.

Setting the heart rate beep

- The pulse rate can be accompanied by a beep when a reading is being taken.
- To turn on the pulse rate beep, navigate to the „Beep" line and hold down the control button to change it to „on". To switch off, proceed in the same way and set to „off".

Adjusting the brightness

- The brightness of the oximeter display can be adjusted.
- To change the brightness, navigate to the „Brightness" line and hold down the control button to set the desired brightness, where „1" is the lowest and „5" the highest.

Factory settings

- If you want to restore factory settings, navigate to the „Restore" line and hold down the control key. „OK" will appear on the line and all lines will be returned to their original values.

5. TAKING A READING

The measurement process itself is very simple. Follow the instructions below.

1. Insert your finger into the finger clamp. Keep your finger still.
2. Press the control button. The oximeter immediately starts measuring. Do not move during the measurement process.
3. After a few seconds, the reading appears on the display.
4. There are 6 different display modes on the device. Use the control button to toggle between them.
5. After removing your finger from the clamp, the oximeter will automatically shut off after 10 seconds.

6. READINGS

SpO2 (oxygen saturation) in%	Classification / Necessary measures
99-94	Normal values
94-90	Lower values: doctor's visit recommended
<90	Critical values: doctor's visit needed

Note: The table does not apply to persons with certain diseases - see Safety instructions (e.g. asthma, heart failure, respiratory diseases) and when staying at altitudes above 1,500 m. If you suffer from any of the diseases mentioned, always consult your doctor, who will assess the readings.

7. TRANSFERRING READINGS TO A MOBILE DEVICE

The TrueLife Oximeter X5 BT allows you to transfer readings to a mobile device and keep them with you at all times. The software automatically saves the measured data and plots the data in the app.

Below you will find instructions on how to transfer the data.

First, download the application, install and register.

1. Download the free application at the link below. The application is intended for mobile devices or other devices running Android or iOS.
2. Install and open the application.
3. Continue by registering in the application. Fill in your email address and password. This information serves as login information.

truelife.eu/56939/app

8. SPECIFICATIONS

Display range	Oxygen saturation	35% to 100%
	Pulse rate	25 bpm to 250 bpm
Resolution	Oxygen saturation	1%
	Pulse rate	1 bpm
Measurement precision	Oxygen saturation	±2% (70% to 100%) No requirement (≤ 69%)
	Pulse rate	±2 bpm
Alert range	Oxygen saturation	Upper limit: 50% to 100% Lower limit: 50% to 100%
	Pulse rate	Upper limit: 25 bpm to 250 bpm Lower limit: 25 bpm to 250 bpm
Alert error	Oxygen saturation	± 1% of the preset value
	Pulse rate	The greater of ±10% of the preset value and ±5 bpm
Battery		High-performance dry battery can be used for about 685 times at normal temperature

EMC INFORMATION – GUIDANCE AND MANUFACTURE'S DECLARATION

1* WARNING: Use of this equipment adjacent to or stacked with other equipment should be avoided because it could result in improper operation. If such use is necessary, this equipment and the other equipment should be observed to verify that they are operating normally.

2* WARNING: Use of accessories, transducers and cables other than those specified or provided by the manufacturer of this equipment could result in increased electromagnetic emissions or decreased electromagnetic immunity of this equipment and result in improper operation.

3* WARNING: Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of the ME equipment, including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result.

Declaration – electromagnetic emission

Emissions test	Compliance
RF emissions CISPR 11	Group 1
RF emissions CISPR 11	Class B
Harmonic emissions IEC 61000-3-2	Not applicable
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Not applicable

Declaration – electromagnetic immunity

Immunity test	IEC 60601 test level	Compliance level
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air
Electrical fast transient/ burst IEC 61000-4-4	± 2 kV for power supply lines ± 1 kV for input/output lines	Not applicable
Surge IEC 61000-4-5	± 0.5kV, ± 1 kV line(s) to lines ± 0.5kV, ± 1 kV, ± 2 kV line(s) to earth	Not applicable
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0 % UT; 0.5 cycle At 0°, 45°, 90°, 135°, 180°, 225°, 270°and 315° 0 % UT; 1 cycle and 70 % UT; 25/30 cycles Single phase: at 0° 0 % UT; 250/300 cycles	Not applicable
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	30 A/m	30 A/m

NOTE: UT is the a.c. mains voltage prior to application of the test level.

Declaration – electromagnetic immunity

Immunity test	IEC 60601 test level	Compliance level
Conducted RF IEC 61000-4-6	3 V 0.15 MHz to 80 MHz 6 V in ISM bands between 0.15 MHz and 80 MHz	Not applicable
Radiated RF IEC 61000-4-3	10V/m 80 MHz to 2.7 GHz	10V/m

Declaration – IMMUNITY to proximity fields from RF wireless communications equipment

Immunity test	IEC60601 test level				Compliance level
	Test frequency	Modulation	Maximum power	Immunity level	
Radiated RF IEC 61000-4-3	385 MHz	**Pulse Modulation: 18Hz	1.8W	27 V/m	27 V/m
	450 MHz	*FM+ 5Hz deviation: 1kHz sine	2 W	28 V/m	28 V/m
	710 MHz 745 MHz 780 MHz	**Pulse Modulation: 217Hz	0.2 W	9 V/m	9 V/m
	810 MHz 870 MHz 930 MHz	**Pulse Modulation: 18Hz	2 W	28 V/m	28 V/m
	1720 MHz 1845 MHz 1970 MHz	**Pulse Modulation: 217Hz	2 W	28 V/m	28 V/m
	2450 MHz	**Pulse Modulation: 217Hz	2 W	28 V/m	28 V/m
	5240 MHz 5500 MHz 5785 MHz	**Pulse Modulation: 217Hz	0.2 W	9 V/m	9 V/m

Note* – As an alternative to FM modulation, 50 % pulse modulation at 18 Hz may be used because while it does not represent actual modulation, it would be worst case.

Note** – The carrier shall be modulated using a 50 % duty cycle square wave signal.

SAFETY WHEN WORKING WITH THE DEVICE:

1. Read the attached safety instructions carefully.
2. For proper handling of the product, read the latest version of the manual, which is available for download at: www.truelife.eu.

Typos and information changes reserved.

Warning: This product is intended for use in households and similar environments. We declare that if the TrueLife Oximeter X5 BT product is used as a professional medical device, a safety inspection (BTK) must be performed regularly at the following intervals: 24 months after the purchase of the product, then every 12 months.

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truelife

EN

Before first use, please read all of the following safety instructions.

General Information
- **Always consult your doctor for specific information about your blood oxygen. Self-diagnosis and treatment using measured results may be dangerous.**
- Read these safety instructions before using the product.
- If a serious adverse event occurs while using the equipment, report it to the manufacturer/importer and to the national regulatory authority.
- For proper handling of the product, read the latest version of the manual, which can be downloaded from the product page at: http://www.truelife.eu.
- Printing errors are reserved.

JUMPER Health App Policies and Notices
- The JUMPER Health app is provided by a third party. Elem6 sro is not responsible for the data provided by you.
- The application is not a medical device. The application is only used to store measured data on a mobile device. The stored measured data are not sent or in any way consulted with doctors.
- The application does not replace a consultation with a specialist, especially a doctor, or any medical care. It is also not used for diagnostic purposes.
- If you have any questions about your medical condition, always consult your doctor.
- Neither the provider of the JUMPER Health application nor the company elem6 sro are providers of health care or medical consultations.

Safety Precautions
- The package contains small parts that can be dangerous for children. Always keep the product out of the reach of children. The bags can cause suffocation if place on the head and the parts they contain are a choking hazard.
- Do not use the device if you are allergic to rubber.
- Do not use the device if the device or the finger used in the device is wet.
- Do not use the device on infants and young children.
- Do not use the device while using a blood pressure monitor.
- Do not use the device during an MRI or CT scan.
- Do not use the device on painted, dirty or manicured fingers.
- Do not use the device on fingers with a large nail thickness, which cannot be inserted into the device without force.
- Do not use the device on fingers with anatomical changes, swelling, scars or burns.
- The device always shows the current measured value, but it cannot be used for continuous monitoring.
- Prolonged use of the device may cause pain in people with circulatory disorders.
- Do not look directly inside the device when taking a reading. Infrared light is harmful to the eyes.
- Keep your hand, finger and body still during the measurement.
- In a person with a heart rhythm disorder, the measured SPO2 and heart rate values may be skewed or the measurement may fail.
- If the device or any component is damaged, do not use it. Contact your dealer or manufacturer.
- For pluggable equipment, the socket-outlet shall be installed near the equipment and shall be easily accessible.

Battery Safety Precautions
- Keep the battery out of the reach of children and animals.
- Use only AAA alkaline batteries for the device.
- Always make sure you insert the batteries with the correct polarity.
- If you are not going to use the device for a long time, remove the batteries.
- Do not throw away worn or damaged batteries in a rubbish bin, fire or heating system, but hand them in at an official hazardous waste collection point.
- Under no circumstances insert any mechanically damaged or swollen batteries into the device. Do not use batteries in this condition at all, there is a risk of explosion.

AC adapter (optional accessory) Safety Precautions
- Use only the AC adapter that is designed specifically for the device. Using any other type of power supply may cause the unit to malfunction or other hazards.
- Do not use the power supply in a highly humid environment. Never touch the power supply with wet hands or while standing in water.
- When powering the device, leave enough space around the charger for air circulation. Do not cover the power supply with paper or other objects that could impair its cooling. Do not use the power supply whilst still in its packaging.
- Connect the power supply to the correct power source. The voltage requirements can be found on the product or on its packaging.
- When disconnecting the AC adapter from the power source, pull the adapter securely. Do not pull on the power cord.
- Do not use the adapter if it is obviously damaged. In case of damage, do not repair it yourself!

Device Maintenance
- Proper maintenance of this device ensures trouble-free operation and reduces the risk of damage.
- Keep the device away from extreme humidity and temperatures, and do not expose the device to direct sunlight or ultraviolet light for extended periods of time.
- Do not drop the device or subject it to strong shocks.
- Do not expose the device to sudden and extreme changes in temperature. This could cause moisture to condense inside the device, which could damage the device. In the event of condensation, allow the device to dry completely before using it again.
- The surface of the display can be easily scratched. Do not touch it with sharp objects.
- Do not clean the device when it is switched on. Use a soft, lint-free cloth to clean the display and the exterior of the device. Do not clean the display with paper towels.
- As a matter of principle, do not attempt to disassemble, repair or modify the device in any way. Disassembly, modification or any attempt at repair will void the warranty and may result in damage to the device or even injury or property damage. Damage to the device may result in inaccurate measurements.
- Do not store or carry flammable liquids, gases, or explosive materials together with the device, its parts, or accessories.
- To prevent theft, do not leave the device and accessories in view in an unattended vehicle.
- Overheating may damage the device.

Manufacturer's Recommendations and Declaration on Electromagnetic Compatibility (EMC)
- Due to the large number of electronic devices, such as personal computers and mobile phones, medical devices may be susceptible to electromagnetic interference from other devices. Electromagnetic interference can cause medical equipment to malfunction and create a potentially hazardous situation.
- In addition, medical devices should not interfere with other devices.
- The EN60601-1-2 standard was used to regulate electromagnetic compatibility (EMC) requirements to prevent hazardous events caused by the product. The standard defines levels of sensitivity to electromagnetic interference and maximum levels of electromagnetic emissions for medical devices.
- This medical device has been manufactured in accordance with EN60601-1-2: 2015. The device should be operated under the following conditions:
- Do not use mobile phones or other devices that generate strong electric and electromagnetic fields near the medical device. This can cause a malfunction and cause a potentially hazardous situation. A minimum distance of 7 m is recommended. Make sure devices are operated at a greater distance.

Further Information
1) For households: The symbol (crossed-out wheelie bin) on the product or in the accompanying documentation means that used electrical and electronic products must not be disposed of with your other household waste. In order to dispose of the product properly, please return it to the designated collection points, where it will be accepted free of charge. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. Contact your local authority or the nearest collection point for further details. Improper disposal of this type of waste may result in penalties in accordance with national regulations. Information for users on disposal of electrical and electronic equipment (corporate and business use): For proper disposal of electrical and electronic equipment, please contact your dealer or supplier for detailed information. Information for users on disposal of old electrical & electronic equipment in countries outside the European Union: The above symbol (crossed-out wheeled bin) is only valid in European Union countries. For proper disposal of electrical and electronic equipment, request detailed information from your authorities or equipment dealer. Everything is expressed by the crossed-out wheeled bin symbol on the product, packaging, or printed materials.

2) Please contact your dealer for warranty repairs. In case of technical problems and questions, please contact your dealer who will inform you of the next steps to be taken. Observe the rules for working with electrical equipment. The user is not authorised to disassemble the equipment or replace any of its components. There is a risk of electric shock if the covers are opened or removed. Incorrect assembly and reconnection of the equipment also exposes you to the risk of electric shock.

The warranty period for the products is 24 months unless otherwise stated. The warranty does not cover damage caused by non-standard use, mechanical damage, exposure to harsh conditions, handling contrary to the manual and normal wear and tear. The warranty period is 24 months for the battery and 6 months for its capacity. For more information on the warranty, please visit www.elem6.com/warranty

Neither the manufacturer, importer nor distributor shall be liable for any damage caused by the installation or improper use of the product.

EU Declaration of Conformity
elem6 s.r.o. hereby declares that the TrueLife Oximeter X5 BT complies with the essential requirements and other relevant provisions of Directive 93/42/ECC and 2014/53/EU. All TrueLife products are intended for sale without restriction in Germany, the Czech Republic, Slovakia, Poland, Hungary and other EU member states. The full Declaration of Conformity can be downloaded from www.truelife.eu/support/doc

Frequency band(s) in which the radio equipment operates: 2.4 GHz ~ 2.48 GHz

Maximum radio frequency power transmitted in the frequency band(s) in which the radio equipment operates: -95 dBm

Printing errors and changes in manual and safety instructions reserved.

Importer:
elem6 s.r.o.,
Braškovská 15, 16100 Prague 6,
Czech Republic

CE 0598 RoHS
Model: JPD-500F

Made in P.R.C.

EC REP MedPath GmbH
Mies-van-der-Rohe-Strasse 8, 80807 Munich, Germany

Shenzhen Jumper Medical Equipment Co., Ltd.
Address: D Building, No. 71, Xintian Road, Fuyong Street, Baoan, Shenzhen, Guangdong, China, 518103

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